



Catalog Number: AP50282HU

Species: Human

Size: 100 µg

Instruction manual

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

---

## Recombinant Human Novel Coronavirus Spike glycoprotein (S) Protein

This package insert must be read in its entirety before using this product.

### **If You Have Problems**

Our expert Technical Support Staff is available to assist you in answering your questions and resolving issues to ensure complete customer satisfaction.

### **Please Contact Us**

**Tel:** (86)-27-65523378

**Fax:** (86)-27-65523378

**Email:** [sales@abebio.com](mailto:sales@abebio.com)  
[service@abebio.com](mailto:service@abebio.com)

**Website:** [www.abebio.com](http://www.abebio.com)

In order to obtain higher efficiency service, please ready to supply the lot number of the kit to us (found on the outside of the box).

### **【BACKGROUND】**

Severe acute respiratory syndrome (SARS) is a viral respiratory illness caused by a coronavirus, called SARS-associated coronavirus (SARS-CoV). Human coronaviruses (HCoVs) were previously only associated with mild diseases. The SARS-CoV genome contains five major open reading frames (ORFs) that encode the replicase polyprotein; the spike (S), envelope (E), and membrane (M) glycoproteins; and the nucleocapsid protein (N).

### **【GENE NAME SYNONYM】**

SARS-CoV-2 Spike glycoprotein Protein; SARS-CoV-2 SgP; Spike glycoprotein protein; novel coronavirus S Protein; novel coronavirus Spike glycoprotein Protein; 2019-nCoV Spike glycoprotein; 2019-nCoV S; 2019nCoV S; 2019-nCoV S Protein; 2019 ncov S Protein; 2019-nCoV Spike glycoprotein protein.

### **【SOURCE】**

Human

### **【HOST】**

*E. coli* 319-541AA.

### **【PROTEIN RESIDUES】**

with N-terminal 6×His-tagged.

### **【PROTEIN SEQUENCES】**

RVQPTESIVRFPNITNLCPFGEVFNATRFASVYAWNRRKRISNCVADYSVLYNSA  
SFSTFKCYGVSPSTKLNLCFTNVYADSFVIRGDEVQRQIAPGQTGKIADYNYKLP  
DDFTGCVIAWNSNNLDSKVGGNYNLYRLFRKSNLKPFERDISTEIYQAGSTP  
CNGVEGFNCYFPLQSYGFQPTNGVGYQPYPYRVVLSFELLHAPATVCGPKKST  
NLVKNKCVNF

### **【PURITY】**

> 90 % as determined by SDS-PAGE.

### **【PREDICTED MOLECULAR MASS】**


**Predicted MW:** 28 kDa

**Observed MW:** 28 kDa

### **【FORMULATION】**

Lyophilized from 10 mM Tris, 1 mM EDTA, 6% Trehalose, pH 8.0

### **【SDS-PAGE】**

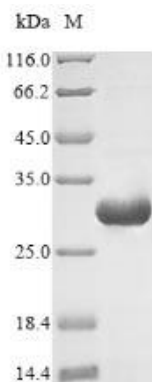

### **【STORAGE】**

Store it under sterile conditions at -20°C to -80°C upon receiving. Recommend to aliquot the protein into smaller quantities for optimal storage.

***\*\*Avoid repeated freeze-thaw cycles.\*\****

### **【STABILITY】**

The recombinant protein is stable for up to 12 months from date of receipt at -80°C.

### **【USAGE】**

**2019 ncov S Protein** - Centrifuge the standard vial at 6000-10000rpm for 30s. We recommend that this vial be briefly centrifuged prior to opening to bring the contents to the bottom. Please reconstitute protein in deionized sterile water to a concentration of 0.1-1.0 mg/mL. We recommend to add 5-50% of glycerol (final concentration) and aliquot for long-term storage at -20°C/-80°C. Our default final concentration of glycerol is 50%. Customers could use it as reference.